R-100135

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R-100135 is a practical Ultrasonic Bath specially designed for a wide range of applications. For a long working life and an optimum performance, we recommend you to read carefully the following instructions:

INSTALLATION

- Carefully unpack the Ultrasonic Bath. Inside you will find the following accessories:
 - 1 plastic basket (F).
 - 1 lid (A).
 - 1 power cord (G).
- Choose a flat and firm surface far from heat or vibration sources.
- Connect the machine to a grounded 230 V, 50/60 Hz power supply.
- Place the basket (F) inside the tank.
- Fill the tank (C) with water or a cleaning solution.
- The machine is ready for operation.

OPERATION

1. Open the lid (A), place the products in the basket (F) in line and then put the basket inside the tank (C).
2. In order to prevent scratches and friction of the products, do not pile them up.
3. Pour water into the stainless steel tank. The water should cover the products but within the *min.* and *max.* marks.
4. Connect the power cord (G) to the body of the unit and the power source. The display will show "180s".
5. 180s is the most common working cycle. You can choose 5 different working cycles (180s→280s→380s→480s→90s→) by pressing the "SET" key.
6. Press the "ON" key to start cleaning. The right green indicator will turn on.
7. During the cleaning, the remaining time is shown on the display.
8. During the cleaning, you can hear the "BUZZ" sound, meaning that the unit is in operation.
9. When time "000" is reached, the unit will stop automatically and the green light will turn off.
10. If another cleaning is needed, just repeat the above steps.

This machine has a protection against overheating. If the right red pilot light turns on, it means that the protection is in operation. If you press the "ON" key, the remaining time will be displayed, but no ultrasonic wave will be produced. The machine needs to stop for 15 minutes. After this time, press the "ON" key and the green pilot light will turn on; the machine can work normally again.

11. Press the "TC" key (J) to start the heating function. The left red pilot light will turn on.
12. If you want to stop the unit during the cleaning, just press the "OFF" key. Press the "ON" key and the unit will continue the cleaning cycle again.
13. To reset the working cycle, just press the "SET" key.
14. When the cleaning is over, unplug the unit, open the lid, and then remove the items. Finally, empty the stainless steel tank and dry it.


CLEANING SOLUTIONS

Cleaning solutions can be classified into two types: those consisting of a cleaning solution mixed with water, and solvents.

The most common cleaning method is with water, as it has an acceptable capacity to dissolve inorganic substances. However, it is extremely limited in dissolving organic substances. For this type of compound, water is a good medium for the dispersion of dirt particles, and a good basis for specific detergents that dissolve organic substances.

In most cases a mixture of water and domestic washing up liquid can be used as a cleaning solution. Nevertheless, when due to the amount or nature of dirt to these cleaning solutions are not enough, specific solvents for ultrasonic baths should be used. **Under no circumstances should inflammable solvents or environmentally dangerous products be used.**

DIRT REMOVING PROCESS

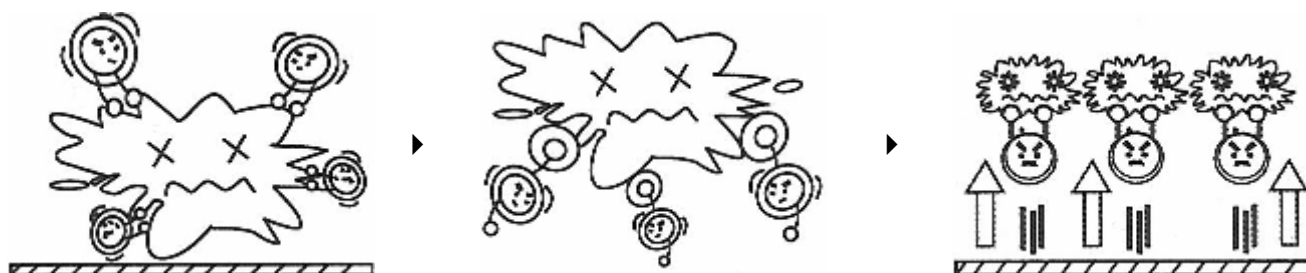The ultrasounds generate myriads of small bubbles.

The ultrasounds make the bubbles implode constantly and impact against the object's surface

The impact breaks up the dirt into many grains and separates them from the object's surface. The grains remain in the water.

MAINTENANCE AND CLEANING

In order to ensure an optimum performance and reliability of the unit, we recommend the following:

- **Avoid excessively wetting the device's keyboard.**
- The unit has several measures to protect the circuit, but if it is operated without water (over 20 seconds) the working life of the machine will be shortened.
- In order to prevent overheating and prolong the machine's working life, it is strongly recommended to stop it for 15 minutes after working for 3 cycles in a row.
- Unplug and empty the unit after using it. Do not keep water in the tank for a long time.
- Dry the tank and unit with a towel. Do not wash the unit completely with water to prevent the risk of shock.
- Do not load the tank excessively; it could reduce the efficacy of the bath and damage the unit.
- Do not put very heavy or large objects in the bottom of the tank.
- Use chemical cleaning products compatible with the objects to be treated. Avoid corrosive solvents.
- When items (especially valuable items) are being cleaned for the first time, carry out a trial run first with just one item. Ask the manufacturer what cleaning solution is recommended.
- During normal operation the bath and the base of the machine can heat up considerably. Temperatures up to 70/80 °C are completely normal, and can be reached after a long and continuous operation with the lid on. If you wish to reduce the temperature, remove the lid and let the machine rest for a few minutes before starting a new cycle.
- After using the unit, clean the tank and the external body with a damp cloth. Always unplug the unit before cleaning it.
- Keep the machine in a cool and dry place.

IMPORTANT- The input voltage of the Ultrasonic Bath transducers can reach hundreds of volts. Therefore, it is very important to observe the following Safety Regulations:

- 1) Do not open the lid when the power cord is connected to the power supply.
- 2) Do not spill water or any other conducting fluid over the unit's housing.
- 3) Always connect the unit to grounded power outlets.

PRECAUTIONS

- ◆ Connect the unit to a grounded 230 V, 50/60 Hz power outlet.
- ◆ Do not use inflammable solvents as cleaning solutions.
- ◆ Keep the unit unplugged before filling it. Do not fill the tank above the *max.* mark. Do not start up the unit when the water level is below the *min.* mark.
- ◆ Do not allow children or unskilled people to handle the unit.
- ◆ Never submerge the unit in any liquid to clean it.
- ◆ Ultrasonic cleaning is not suitable for organic pearls.
- ◆ Do not open the unit while it is connected to the power supply.

- ◆ The machine lid may be hot after a cycle. Proceed with caution when removing it.
- ◆ If you are not sure of the effects of the cleaning solution on a given type of item, carry out a previous test with a single item.
- ◆ Do not attempt to remove any screw or panel. The user does not need to manipulate anything inside.
- ◆ To prevent electric shock, never touch the power cord or socket with wet hands.
- ◆ Never operate the unit if the plug or the power cord is damaged. If the unit does not work properly, is damaged or has fallen into the water, refer it to a technical service for examination and repair.
- ◆ Do not touch the unit if it has fallen into the water. Unplug it immediately.
- ◆ Use the unit properly according to the user's guide. Do not use accessories not recommended by the manufacturer.
- ◆ The unit should never be left unattended when it is plugged.
- ◆ Do not operate the unit next to heat sources such as radiators.
- ◆ Install the unit on a dry and flat surface.

TECHNICAL SPECIFICATIONS

Voltage:	230 V, 50/60 Hz
Power:	170 W
Frequency:	42 KHz
Tank material:	Stainless steel
Tank dimensions:	264 x 164 x 80 mm
Outside dimensions:	340 x 255 x 225 mm
Tank capacity:	2500 ml <i>min:</i> 600 ml / <i>max:</i> 2100 ml
Timer settings:	180s→280s→380s→480s→90s
Weight:	2.5 kg

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***TALLERES MESTRAITUA S.L.** is a company manufacturing tools and appliances for the dental prothesis sector, from the business address shown at the bottom of the page.*

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***I HEREBY STATE:** That the **ULTRASONIC CLEANING BATH R-100135**, with the following characteristics:*

Tensión de alimentación / <i>Voltage:</i>	230 V, 50/60 Hz
Potencia / <i>Power:</i>	170 W
Frecuencia / <i>Frequency:</i>	42 kHz
Material de la cuba / <i>Tank material:</i>	Acero inoxidable / <i>Stainless steel</i>
Dimensiones de la cuba / <i>Tank dimensions:</i>	264 x 164 x 80 mm
Dimensiones exteriores / <i>Outside dimensions:</i>	340 x 255 x 225 mm
Capacidad de la cuba / <i>Tank capacity:</i>	2100 ml
Posiciones del temporizador / <i>Timer settings:</i>	180s→280s→380s→480s→90s
Peso / <i>Weight:</i>	2.5 kg

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Meets the low voltage requirements set out in Directive 2006/95/EC and in Directive 2004/108/EC concerning electromagnetic compatibility.

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In witness whereof, the company issues this declaration of conformity in Sondika, on March 12th 2010.

Firma / *Signed:*

Ignacio Mestraitua,
Director Gerente / *Managing Director*

MESTRA®

TALLERES MESTRAITUA S.L.

Txori-Erri Etorbidea, 60

Tfno. + 34 944530388 - Fax + 34 944711725

E-mail: mestra@mestra.es - www.mestra.es

48150 SONDIKA - BILBAO - ESPAÑA